Product Datasheet

NFkB-p65(Ab-276) Antibody

Catalog No: #21011



Package Size: #21011-1 50ul #21011-2 100ul #21011-4 25ul


Overview

Product Name	NFkB-p65(Ab-276) Antibody
Host Species	Rabbit
Clonality	Polyclonal
Applications	WB IHC IF
Species Reactivity	Hu Ms Rt
Immunogen Type	Peptide-KLH
Target Name	NFkB-p65
Alternative Names	NFKB3; RELA; TF65; Transcription factor p65; p65


Application Details

Predicted MW: 65kd	
Western blotting: 1:500~1:1000	
Immunohistochemistry: 1:50~1:100	
Immunofluorescence: 1:100~1:200	


Images

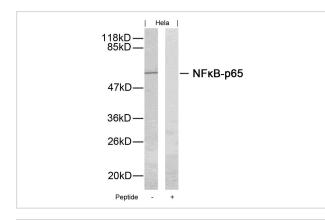

Western blot analysis of extracts from Hela cells using NFkB-p65(Ab-276) Antibody #21011 and the same antibody preincubated with blocking peptide.

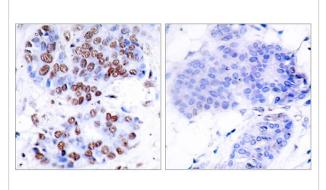

Immunohistochemical analysis of paraffin-embedded human breast carcinoma tissue using NFkB-p65(Ab-276) Antibody #21011(left) or the same antibody preincubated with blocking peptide(right).

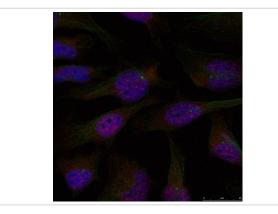

Immunofluorescence staining of methanol-fixed Hela cells using NFkB-p65(Ab-276) Antibody #21011.

Descriptions	
Immunogen	Peptide sequence around aa.274~278 (R-P-S-D-R) derived from Human NFkB-p65.
Specificity	The antibody detects endogenous level of total NFkB-p65 protein.
Purifiction	Antibodies were produced by immunizing rabbits with synthetic peptide and KLH conjugates. Antibodies were
	purified by affinity-chromatography using epitope-specific peptide.
Formulation	Supplied at 1.0mg/mL in phosphate buffered saline (without Mg2+ and Ca2+), pH 7.4, 150mM NaCl, 0.02%
	sodium azide and 50% glycerol.
Storage	Store at -20°C for long term preservation (recommended). Store at 4°C for short term use.
Accession NO.	Swiss-Prot: Q04206NCBI Protein: NP_001138610.1

Related Information

NF-kappa-B is a pleiotropic transcription factor which is present in almost all cell types and is involved in many biological processed such as inflammation, immunity, differentiation, cell growth, tumorigenesis and apoptosis. NF-kappa-B is a homo- or heterodimeric complex formed by the Rel-like domain-containing proteins RELA/p65, RELB, NFKB1/p105, NFKB1/p50, REL and NFKB2/p52 and the heterodimeric p65-p50 complex appears to be most abundant one. The dimers bind at kappa-B sites in the DNA of their target genes and the individual dimers have distinct preferences for different kappa-B sites that they can bind with distinguishable affinity and specificity. Different dimer combinations act as transcriptional activators or repressors, respectively. NF-kappa-B is controlled by various mechanisms of post-translational modification and subcellular compartmentalization as well as by interactions with other cofactors or corepressors. NF-kappa-B complexes are held in the cytoplasm in an inactive state complexed with members of the NF-kappa-B inhibitor (I-kappa-B) family. In a conventional activator pathway, I-kappa-B is phosphorylated by I-kappa-B kinases (IKKs) in response to different activators, subsequently degraded thus liberating the active NF-kappa-B p65-p65 complex appears to be involved in invasin-mediated activation of IL-8 expression. The inhibitory effect of I-kappa-B upon NF-kappa-B the cytoplasm is exerted primarily through the interaction with p65. p65 shows a weak DNA-binding site which could contribute directly to DNA binding in the NF-kappa-B complex

Baeuerle P A, et al. (1994) Annu Rev Immunol. 12:141-179. Baeuerle P A, et al. (1996) Cell 87:13-20. Haskill S, et al. (1991) Cell 65:1281-1289. Thompson J E, et al. (1995) Cell 80: 573-582.

Published Papers

Gerben M. van Buul, Wendy L.M. Koevoet, Nicole Kops el at., Platelet-Rich Plasma Releasate Inhibits Inflammatory Processes in Osteoarthritic Chondrocytes., Am J Sports Med., 39: 2362-2370(2011)

PMID:21856929

S. Clockaerts, Y.M. Bastiaansen-Jenniskens, C. Feijt el at., Peroxisome proliferator activated receptor alpha activation decreases inflammatory and destructive responses in osteoarthritic cartilage., Osteoarthritis and Cartilage, 19(7):895-902(2011)

PMID:21458577

Note: This product is for in vitro research use only and is not intended for use in humans or animals.